


Angiographic Catheter


Precision, Performance, Perfection.

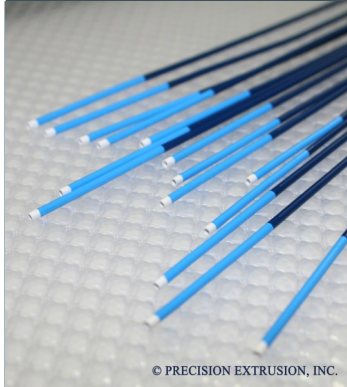

Optimal use of unique resins, radiopaque fillers, and braid wire suitable for medial use provides:

- ✧ Superior pushability
- ✧ Excellent torque control
- ✧ Ease of curve placement
- ✧ Outstanding curve retention
- ✧ Kink resistance
- ✧ Injection stability - Optimal back-out resistance
- ✧ Smooth lumens
- ✧ High pressure
- ✧ Excellent catheter visualization

- ✧ **Assembly design:** allows for easy curve setting – tapering and perfusion port forming when required
- ✧ **Standard French sizes:** providing ease of use with a wide variety of sheath introducers
- ✧ **Encapsulated stainless steel wire braid:** providing kink resistance, super torque and increased burst strength during power injections
- ✧ **Manufacturing process:** providing reliability, super smooth surface and ease of use with introducer
- ✧ **Catheter Curve Setting Technology:** comes from years of experience and multiple doctors feedback
 - Proprietary processing
 - Detailed drawings and curve forms for 35 common curves for 4, 5, and 6F
 - Perfusion ports as necessary

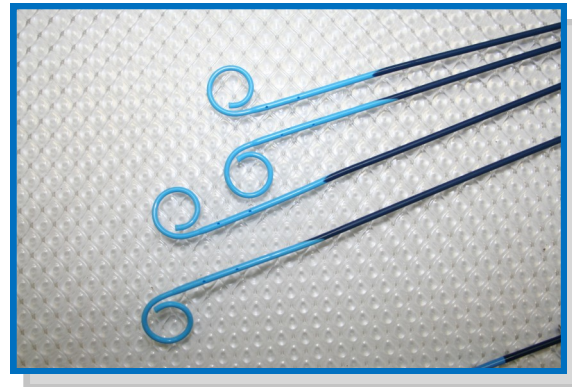

- ✧ **Exceptional flow rates:** ease of injection due to large lumen and selected materials
- ✧ **Variable resins:** utilized allowing for variable stiffness, great pushability, injection

Contact us:

12 Glens Falls Technical Park, Glens Falls, New York 12801 USA

E-mail: sales@precisionextrusion.com **Web Site:** www.precisionextrusion.com **Phone:** +1 (877)-792-1199

Angiographic Catheter

Precision, Performance, Perfection.

Functional lumen sizes and large lumens provide:

- ✧ Exceptional flow rates
- ✧ Ease of hand injection

PEI Angiography Catheter Specifications - Select

Size	4FR Select	5FR Select	6FR Select
ID Shaft (Braided Section) (+/-0 .001")	0.0420"	0.0465"	0.0565"
ID Soft Ext/Tip Select I & II (+/-0 .001")	0.0420"	0.0465"	0.0565"
OD Shaft/Soft Ext/Tip - Select I & II (+/-0 .0010"-0.0015")	0.0540"	0.0660"	0.0790"
Length Soft Ext - Select I (+/- 2mm)	28mm	35mm	35mm
Length Soft Ext - Select II (+/- 2mm)	75mm	80mm	80mm
Length Tip - Select I & II (+/- 0.5mm)	2mm	2mm	2mm
Welded Shaft Length (+/-4cm)	104cm	104cm	104cm

PEI Angiography Catheter Specifications - Pigtail

Size	4FR Pigtail	5FR Pigtail	6FR Pigtail
ID Shaft (Braided Section) (+/-0.001")	0.0420"	0.0465"	0.0540"
ID Soft Ext non tapered section +/-0 .0010"	0.0420"	0.0465"	0.0520"
OD Soft Ext non tapered +/-0 .0010"-0.0015" section	0.0540"	0.0660"	0.0790"
ID Soft Ext tapered (+/-0.002") section	not tapered	0.0430"	0.0440"
OD Soft Ext Tapered (+/-0.002") Section	not tapered	0.0600"	0.0630"
Length Tapered Portion (+/-2mm)	not tapered	45mm	48mm
Length Overall Pigtail Soft Ext (+/-2mm)	75mm	80mm	80mm
Welded Shaft Length (+1/-1.5cm)	111.5cm	111.5cm	111.5cm

Contact us:

12 Glens Falls Technical Park, Glens Falls, New York 12801 USA

E-mail: sales@precisionextrusion.com **Web Site:** www.precisionextrusion.com **Phone:** +1 (877)-792-1199